

# Active 2

VERSATILITY AND COMFORT



The Softform Premier Active 2 retains properties as a high risk reactive support surface however should a patient require stepping up to an Active support surface, then a discrete pump can be fitted transforming the mattress into a dynamic surface and delivering additional levels of pressure relief.



Yes, you can.®

Making Life's Experiences Possible





# Softform Premier Active 2

## Unique Features

## Benefits

<b>Air Insert contained within standard Softform Premier mattress foam construction</b>	Using the pump allows stepping up and down between Reactive and Active alternating systems with minimum patient disruption
<b>Inflation of air cells alternates on a 10 minute cycle</b>	Promotes tissue recovery and pressure distribution
<b>New air cell design</b>	For the users safety and peace of mind
<b>Air system placed beneath foam insert</b>	Greatly reduces risk of sharps puncture
<b>Mattress can stay on the bed for a wide range of patients.</b>	Reduces moving and handling issues; minimizes additional storage requirements.
<b>High risk static mattress, delivering exceptional pressure reduction and patient comfort</b>	In the event of a power failure, the patient is still resting on a high risk pressure reducing surface, minimizing the potential of pressure damage occurring until power is restored
<b>Mattress does not require turning Reduces manual handling issues</b>	Reduces manual handling issues
<b>Two- way stretch PU waterproof cover</b>	Provides excellent recovery and durability, and helps reduce shear and friction forces
<b>High frequency welded seams</b>	Reduces the risk of fluid ingress
<b>Full length flap conceals the zip</b>	Reduces the possible fluid ingress through the zip
<b>Toughened PU coated base</b>	Prolongs the longevity of the mattress
<b>Single piece of foam "U" core</b>	Facilitates patient transfer by providing side wall support
<b>Water penetration resistance</b>	35 kPa minimum (BS 3424: part 26). 35 kPa minimum after 20,000 cycles Martindale Abrasion Test EN ISO 5470-2

## Tech Specs




### WATERFLOW SCALE


Waterflow Scale  
25




User weight  
247.6kg

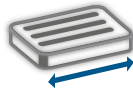

Length  
197cm

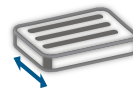

Width  
88cm

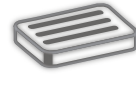

Depth  
15.2cm

### WARRANTY

Warranty  
4yr cover  
8yr foam  
2yr pump

### FOR FURTHER INFORMATION:

**New Zealand:** Freephone 0800 468 222, email [sales@invacare.co.nz](mailto:sales@invacare.co.nz) or visit [www.invacare.co.nz](http://www.invacare.co.nz)

**Australia:** Freephone 0800 460 460, email [sales@invacare.com.au](mailto:sales@invacare.com.au) or visit [www.invacare.com.au](http://www.invacare.com.au)

**Asia:** Phone +61 2 8839 5330, email [exportsales@invacare.com.au](mailto:exportsales@invacare.com.au) or visit [www.invacare.com](http://www.invacare.com)

**Thailand:** Phone +66 0 2108 8555 or visit [www.invacare.com](http://www.invacare.com)



Yes, you can.®