

**permobil**

ROHO SEATING + POSITIONING

SELECT SERIES  
CUSHIONS



**EXTRAORDINARY.  
JUST LIKE YOU.**

# SUPERIOR SKIN & SOFT TISSUE PROTECTION HAS NEVER FELT BETTER

## The SELECT Series

The SELECT Series is one of ROHO's many popular product lines, with the revolutionary ISOFLO® Memory Control which offers advanced skin and soft tissue protection as well as positioning capabilities. The ISOFLO provides each individual an adjustable, quick and simple stability and positioning solution while seated.

### BENEFITS

#### Skin & Soft Tissue Protection

Based on over 45 years of research and individual experience, ROHO® DRY FLOATATION® Technology provides unparalleled performance in facilitating blood flow by fitting and matching the shape of an individual.

Select Series cushions are one of our most popular product lines because of the revolutionary ISOFLO Memory Control, which offers shape fitting capabilities with the push of a knob. The ISOFLO Memory Control offers a seated individual an adjustable, quick and simple positioning and stability solution.

"I love my QUADTRO SELECT Cushion and think it is pretty cool! It helps me sit in my wheelchair and protects my skin."

Abel Rose – Waxahachie, TX



### FEATURES

#### The perfect sitting position is only a push away

The ISOFLO Memory Control, located at the front of all SELECT Series cushions, allows air to flow through the cushion in two ways:



With the ISOFLO unlocked, air can flow freely throughout the entire cushion, similar to a ROHO Single Compartment Cushion.



With the ISOFLO locked, air is contained in four separate quadrants offering stability and positioning capabilities.

## PRODUCTS



### QUADTRO SELECT® HIGH PROFILE® Cushion

Approximate Height: 4.25 inches, 10.5 cm

Approximate Immersion: 3.75 inches, 9 cm

Accommodation of pelvic asymmetry, up to approximately 3 inches, 7.5 cm



### QUADTRO SELECT LOW PROFILE® Cushion

Approximate Height: 2.25 inches, 5.5 cm

Approximate Immersion: 1.75 inches, 4 cm

Accommodation of pelvic asymmetry, approximately up to 1 inch, 2.5 cm




### CONTOUR SELECT® Cushion

Approximate Height (Multiple cell height): 2.25 – 4.25 inches, 5.5 – 10.5 cm

Approximate Immersion: 1.75 inches, 4 cm

Accommodation of pelvic asymmetry, up to approximately 1 inch, 2.5cm

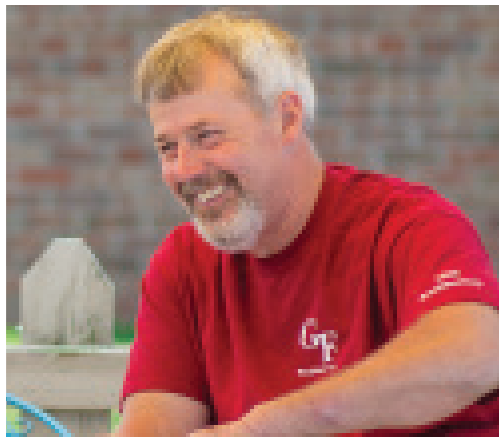

## WHY ROHO?

### Protection

We are the worldwide leader of seating solutions to prevent and treat pressure injuries including deep tissue injury while providing skin and soft tissue protection.

### Proof

ROHO has been in over 90 scientific and clinical studies to reinforce and verify that a ROHO DRY FLOATATION® cushion heals, treats and prevents pressure injuries including deep tissue injury.

### Adjustability

ROHO DRY FLOATATION® cushions, unlike foam products, provide the ability to constantly adjust and conform to an individual's body, now, and as the individual changes throughout the months and years.

### Adaptability

ROHO DRY FLOATATION® cushions support the shifts and movements that individuals make throughout each day.

### Durability

The high quality of ROHO's proprietary neoprene and our ISO certified manufacturing processes instill confidence that the cushion will continue to perform over years.



## ROHO QUADTRO SELECT® CUSHIONS

### APPLICATIONS

Based on clinical, scientific or engineering evidence, this product may be suitable for individuals who:

- have a history of skin/soft tissue breakdown and/or any category/stage pressure injury (including deep tissue injury), skin grafts, and/or myocutaneous flap surgery
- currently have skin/soft tissue breakdown and/or any category/stage pressure injury (including deep tissue injury)
- are at risk for skin/soft tissue breakdown
- have normal or impaired sensation
- require increased stability
- require correction of pelvic asymmetry (medial-lateral and/or anterior-posterior)
- require different positioning of the left and right femur

### SPECIFICATIONS

**INCLUDES:** QUADTRO SELECT HIGH PROFILE, MID PROFILE or LOW PROFILE Cushion, two-way stretch cover, hand inflation pump, repair kit, and operation manual.

**CONSTRUCTION:** Flame-resistant, neoprene rubber (made without latex).

**SIZES\*:** Fits chair width

HIGH PROFILE - 11" to 25"/28 cm to 63.5 cm

MID PROFILE - 13" to 19"/33 cm to 48 cm

LOW PROFILE - 10" to 24"/25.5 cm to 61 cm

\*Customization and custom sizes available

**WEIGHT:** Varies by size, approximately 1.6 kg

**WEIGHT LIMIT:** Unlimited (Cushion must be properly sized to the individual)

**CUSHION WARRANTY:** 36-Month Limited Warranty

**COVER WARRANTY:** 6-Month Limited Warranty

## ROHO CONTOUR SELECT® CUSHION

### APPLICATIONS

Based on clinical, scientific or engineering evidence, this product may be suitable for both the individuals listed under the QUADTRO SELECT Cushions and for individuals who:

- benefit from a pre-contoured seating surface
- require lower extremity alignment
- require a greater amount of pelvic control – the cushion's pre-contoured surface allows for aggressive positioning
- require accommodation of pelvic asymmetry, up to approximately 1 in./2.5cm

### SPECIFICATIONS

**INCLUDES:** CONTOUR SELECT Cushion, two-way stretch cover, hand inflation pump, repair kit, and operation manual.

**CONSTRUCTION:** Flame-resistant, neoprene rubber (made without latex)

**SIZES\*:** Fits chair width – 13" to 24"/33 cm to 61 cm

\*Customization and custom sizes available

**WEIGHT:** Varies by size, approximately 1.4kg

**WEIGHT LIMIT:** Unlimited (Cushion must be properly sized to the individual)

**CUSHION WARRANTY:** 36-Month Limited Warranty

**COVER WARRANTY:** 6-Month Limited Warranty

## WHY PERMOBIL?



### Expertise

Permobil, TiLite, and ROHO®. Three brands. One Mission. For over 45 years, Permobil has been a leader in the complex rehabilitation power wheelchair industry. ROHO became incorporated in 1973, and TiLite has been around for almost three decades, with custom manual wheelchairs being fabricated since the 1990s. Together there is over 90 years of business experience in complex seating and mobility products. And together, we are better able to innovate and create healthcare solutions.

By our commitment to using technology to meet the user's needs, the Permobil Group continues to lead the industry with a diverse line of innovative, complex rehabilitation solutions. This includes the user-focused process of custom fitting and configuration of TiLite manual wheelchairs, the ground breaking anterior tilt technology of the Permobil F-Series, and the decades of research demonstrating the effectiveness of ROHO DRY FLOTATION® Technology. It's not just the "power of mobility" anymore, it's a broad commitment to healthcare solutions.

©2016 ROHO, Inc. The following are trademarks and registered trademarks of ROHO, Inc:

CONTOUR SELECT®, DRY FLOTATION®, HIGH PROFILE®, ISOFLO®,

LOW PROFILE®, MID PROFILE™, QUADTRO SELECT®, ROHO® and shape fitting technology®.

We recommend consulting with a certified rehab professional when ordering and fitting ROHO products.

**permobil**

ROHO SEATING + POSITIONING



The Quality Management System of ROHO, Inc. is certified to ISO 9001:2008 and ISO 13485:2003.